
TRUESONIC

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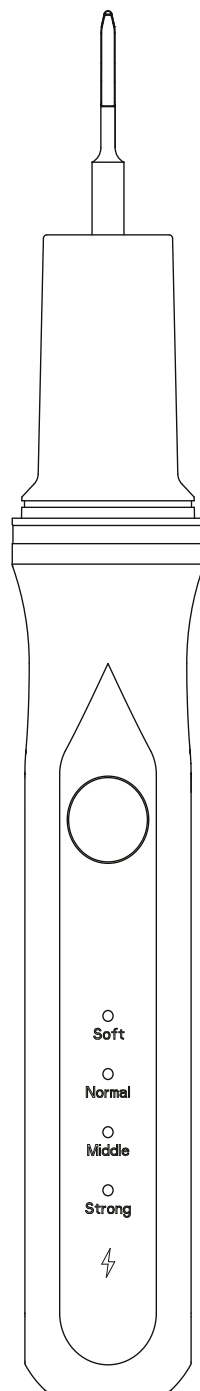Product introduction

Health starts with your teeth. Thank you for choosing our Truesonic Intelligent Ultrasonic Toothbrush. It is specially designed for home use. Truesonic has a stable performance with a variety of functions. It effectively removes plaque, tartar, stains from teeth and prevents gum disease by removing hard tartar. This multi-functional plaque remover uses advanced vibration to remove tartar and pigmentation from your teeth. Operation is simple, comfortable, fast and effective.

How to use it safely

- Before using the product, please read this manual carefully AND KEEP IT FOR FUTURE USE.
- Always place the tip against the surface of the tooth when using. Avoid unnecessary contact of the tip with soft tissues such as gums, mucous membranes or skin.
- You will only feel the vibration when it comes into contact with the tooth. The cleaner will automatically switch off when it comes into contact with the gum.
- The tip must not touch the ceramic denture, filling or any tooth filling during cleaning. The tip could be damaged. Do not use the device by any method other than that described in this manual.
- The metal head of the tooth cleaner is replaceable.
- Keep the toothbrush well away from the eyes.
- In case of malfunction, please stop using.
- This product is not a toy. Do not let your children use it.
- Forbidden to patients with pacemaker.
- Patients with heart disease, pregnant women and children should exercise caution.
- Do not expose the product to direct sunlight, it should be stored 0° to 40°C indoors.
- Do not disassemble the product yourself.
- Do not use the product while it is being charged.
- In case of quality problems with the product, contact the customer service of the store where you purchased the product.

Product principle

This cleaner works on the principle that converts electrical energy into ultrasonic vibration energy. Press the power button to switch on the device, adjust the power using the mode buttons and start cleaning the desired areas on the teeth. The high vibration is used to crush tartar, remove stubborn stains and loose tartar.


Matters needing attention

1. If you have periodontal disease, consult your dentist before use.
2. Bleeding may occur even if the gums are healthy.
3. If gum bleeding persists, which may be caused by gum defects such as periodontal disease, discontinue use and consult your dentist.
4. Do not leave the working tip on the surface of the teeth for a long time during the brushing process.
5. No sound of the tip is heard when the device is started. You will only hear the ultrasonic sound when it comes into contact with the tooth.
6. Do not soak in water or liquid.
7. Keep the machine clean before and after use.
8. The working tip must be tightened when in use.
9. Do not bend or grind the working point.
10. Do not knock or scratch the handle.
11. Dry the machine after use and store it in a cool, dry environment.
12. Do not pull on the power cord while using the machine to avoid separating the handle from the machine.

Cleaning and maintenance

1. Do not insert or drop the product into or into other liquids.
2. Handle the equipment with care and store it in a dry place.
3. Do not mix with toxic, corrosive, flammable or explosive materials.
4. Rinse your mouth thoroughly after each cleaning.
5. Wash the metal head under clean running water.
6. 6. Wipe the body of the toothbrush.


Product description

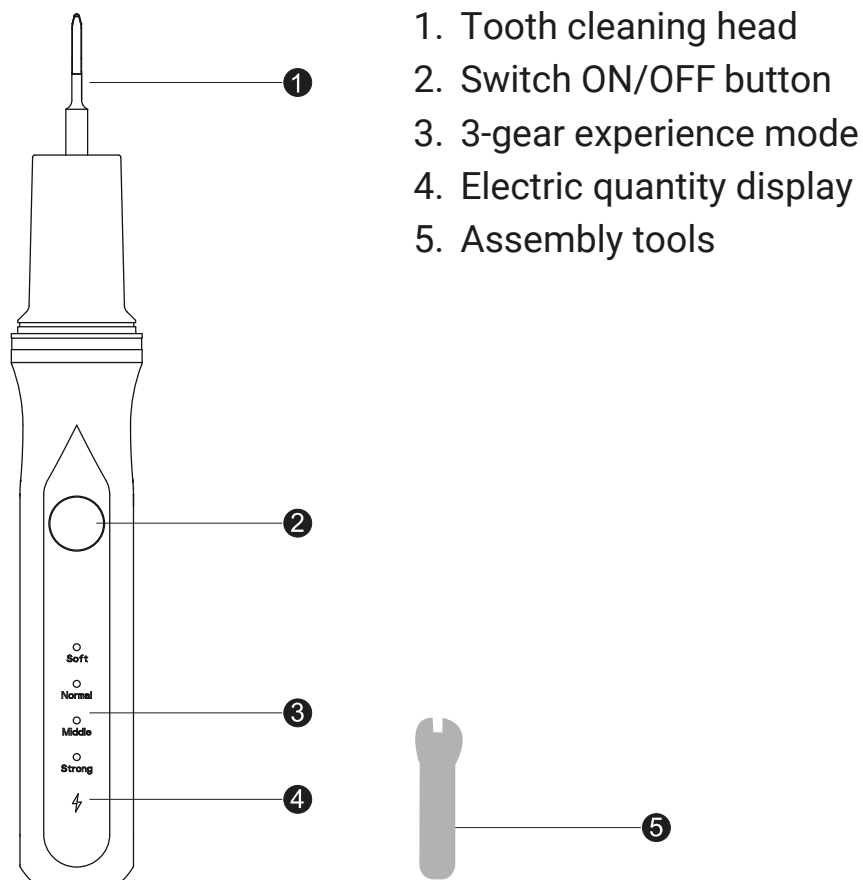
How to operate

1. Make sure the cleaner is fully charged before use. It takes approximately 2 hours to fully charge.
2. If this is the first time you have charged the toothbrush or if you have not used it for a long time (more than 3 months), it may take longer to fully charge.
3. Select a suitable tip and install it clockwise. DO NOT use strong rotation to avoid damage.
4. Press the button to start the unit ON and turn on the SOFT mode for normal cleaning. The start mode indicator light and LED light on the front end of the handle will illuminate, indicating the machine is running.
5. Place the side of the tip end on the tooth so that the tip movement is parallel to the tooth surface.
6. However, difficult to remove tartar may require a higher power setting in COMFORT and STRONG mode.
7. Be sure to maintain tip contact with the tooth surface during cleaning. Keep moving the tip slowly back and forth and let the machine do the work. Use short and long strokes to clean the entire surface of each tooth.
8. Never linger in one place for more than 5s to avoid heating up the treatment area.
9. Caution: The device is susceptible to damage to internal and external components caused by dropping or improper handling.
11. Use the appliance with a moistened mouth.
12. If your mouth dries out during use, rinse it out.
13. When using on the inside of a tooth, use a dental mirror.
14. Low voltage detection. The alarm flashes when the battery voltage is below 3V.
15. After use, long press the OFF button and turn off the power. The device has a preset auto shut-off function after 10 minutes to prevent overheating protection.

Introduction to work tips

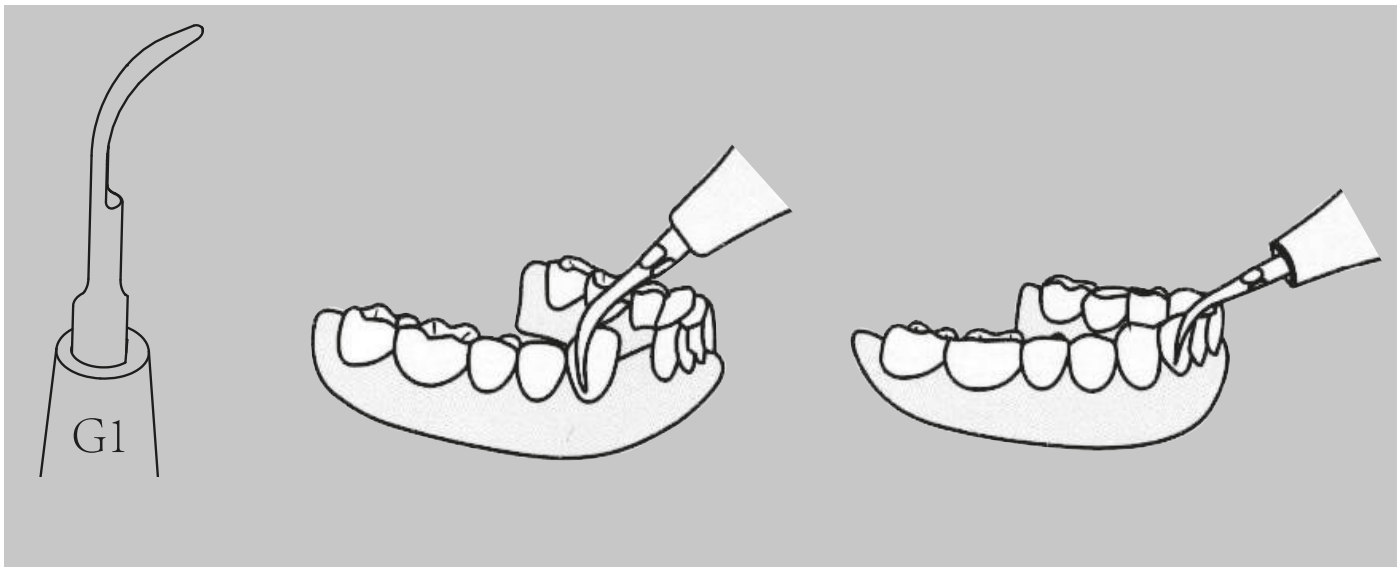Tip 1 is thin, used for tartar removal and interdental cleaning.

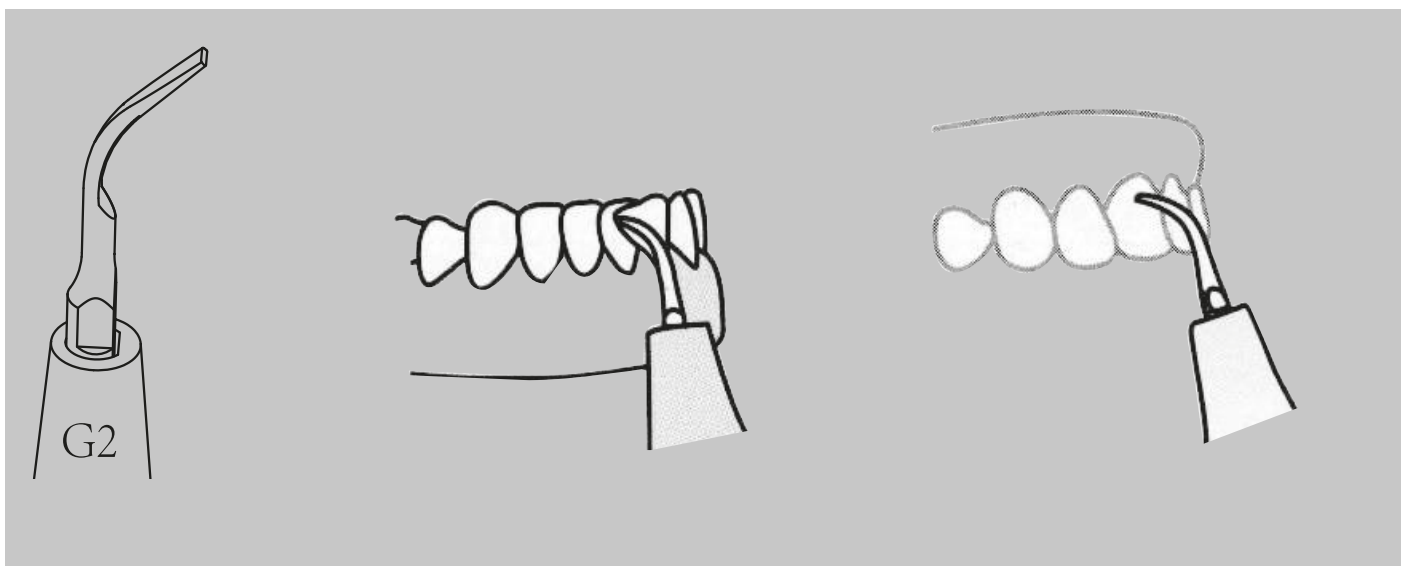
Tip 2 is wider, used to remove large tartar and pigmentation on the teeth.

- Disinfect the working tip with hot water at high temperature before use.
- To prevent cross-infection, do not share the working tip with other people.
- The frequency of the instrument is different, please choose the correct working mode.
- Please use the working tip correctly as above.



How to use work tips

The working tip gently touches the teeth vertically and moves up and down as well as left and right. Adjust the time of use of each part according to the condition of the teeth.


⚠ Please do not use the surface of the work tip without tooth dirt.

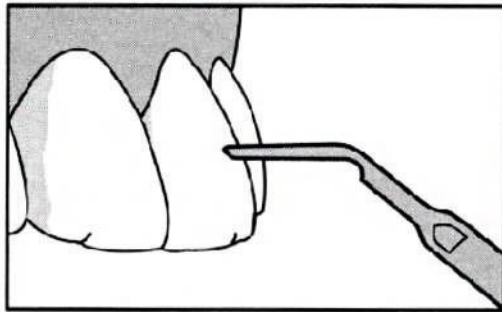
Tooth surface dirt

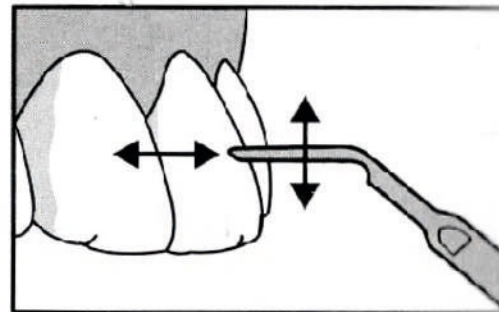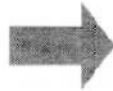
Boundary between gum and tooth

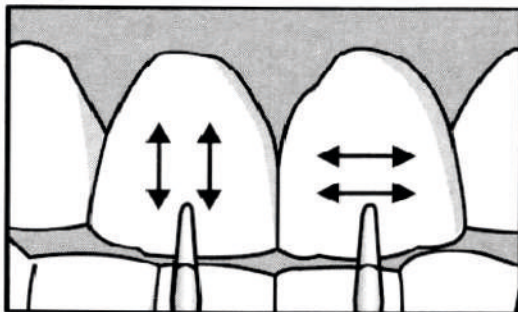
Interdental gap


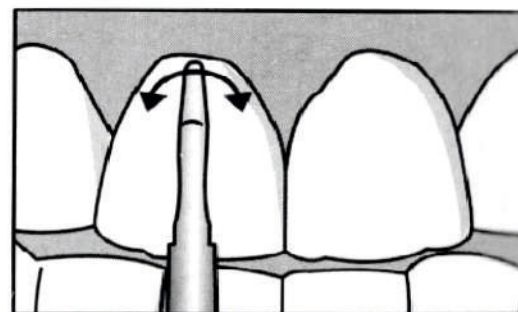
Dirt from uneven teeth

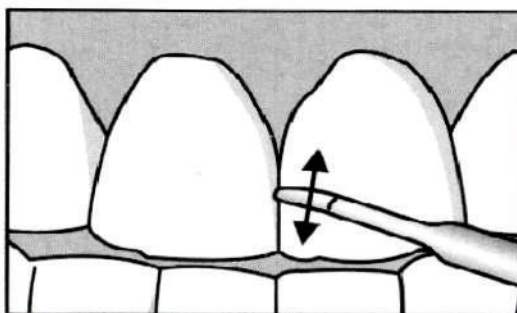Dirt on the back of teeth

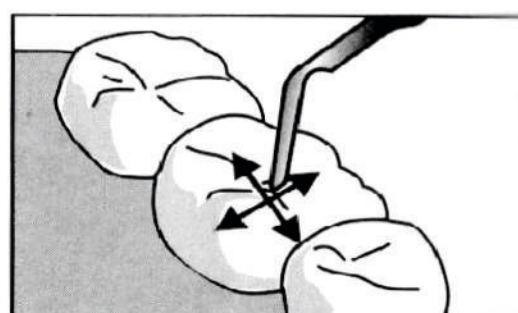Please use dental glass when treating the back of teeth

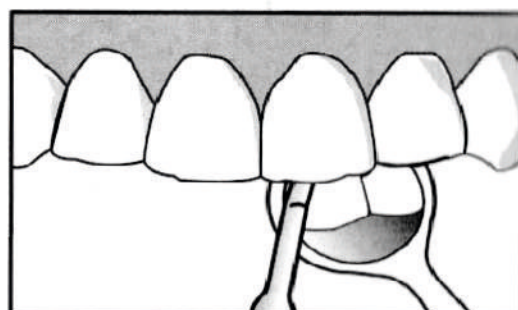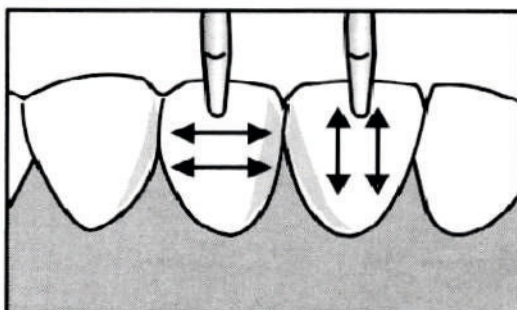
How to change the work tip


The orientation of the work tip may be different but as long as it is firmly mounted on it, it can be used in any direction.


If the connection is loose, vibrations cannot be transmitted properly. Use the installation wrench and make sure the working tip is connected firmly. Do not use other tools to disassemble. You may cause malfunction or damage.

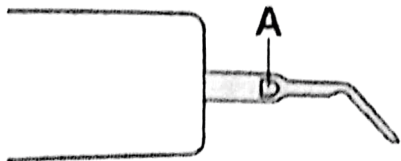
Confirm position A

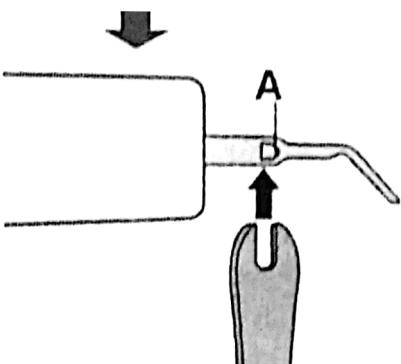Insert the wrench in position A

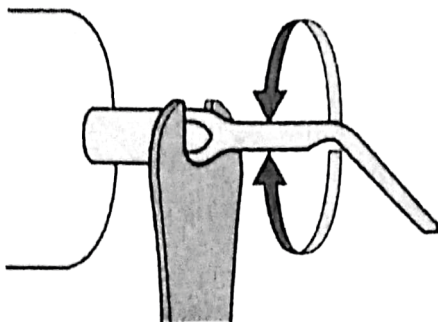Rotate the head counter clockwise with a wrench to separate the working tip and complete the disassembly.

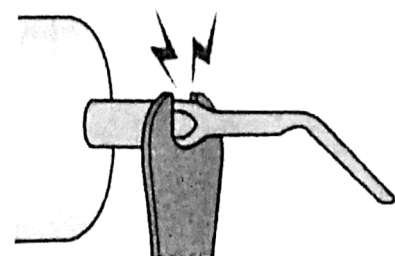
Preventive measure

Avoid over tight assembly, resulting in damage to the interface.

How to maintain

Please do not use diluents, benzene and solvents, chemicals, scrubbing powder.

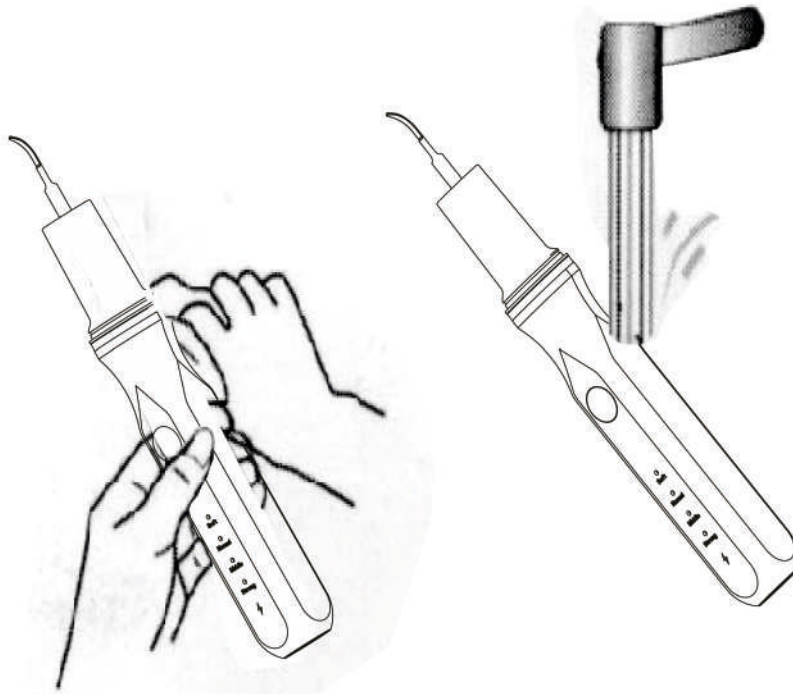The main body can be directly washed with water, and the dirt on the main body can be wiped with a rag.

Notice

- Please consult a licensed dentist or dental hygienist before using our appliances. We provide a product, not a treatment plan, and are therefore not responsible for any complications. If you choose to use it without professional supervision, you use it at your own risk. When you contact us and if we give you any advice, do not take it for granted as we are not responsible for the advice we have given you.
- This product is not a substitute for medical equipment and has no therapeutic effect.
- The product is not suitable for children under 18 years of age.